Double Sheathed Dispensing Needle

Suitable for the reconstitution and preparation of freeze-dried powder drugs.

Version 1.1_2024-11-06



Protective Cap

External Needle

Aspirate the liquid from the ampoule

Internal Needle

Puncture the rubber stopper of vial to reconstitute and aspirate the drug liquid



Incline-hole Needle

Aspirate less drug liquid with high viscosity




Side-hole Needle

Aspirate more drug liquid with low viscosity

Product Features

- Uniquely patented double sheathed needle design
- Silicification of needle tip can decrease the puncture resistance and reduce operation intensity of medical staff
- Made of high-strength materials, with low risk of puncture fracture
- Made of polymer material, it can decrease the risk of accidental puncture and reduce the puncture chips

Procedures

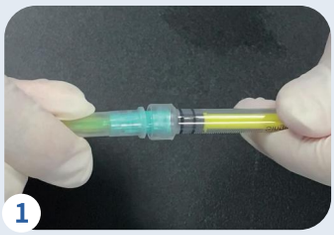

1 Tighten the double sheathed dispensing needle to the syringe



2 Insert the yellow external needle into the ampoule until the end and aspirate the diluent




3 Twist and pull out the yellow external needle



4 Inject the diluent into the vial and gently shake the vial



5 Invert the vial and slowly aspirate the solution

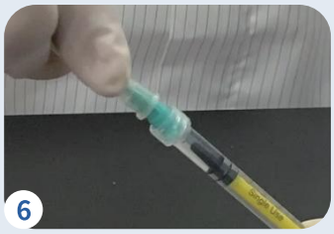

6 After covering the protective sheath, screw off the internal needle and discard it

* This product can only be operated by qualified medical personnel, and used in medical institutions.

Nextech-X

Hangzhou Qiantang Longyue Biotechnology Co., Ltd.
Rooms 302 and 104, Building 12, Block 1, No.619 Wangmei Road, Linping Street, Linping District, Hangzhou City
Website: www.nextech-x.com/en
Email: Sales1@Nextech-x.com
Tel: 0571-89150121, 18017582758