Nasal Atomized Device

This product atomizes the liquid drug for rapid and painless clinical nasal atomized administration.

Version 2.0_2025-06-12





Product Features

- *Needle-free and painless*
- *Rapid absorption*
- *Flexible administration*
- *Patient and paramedic friendly*
- *Micron atomization*



Syringe



Aspirating Needle



Nozzle




Nozzle (soft material)

Nozzle angle can be adjusted among 180°


* This product can only be operated by qualified medical personnel, and used in medical institutions.

Nextech-X


Procedures

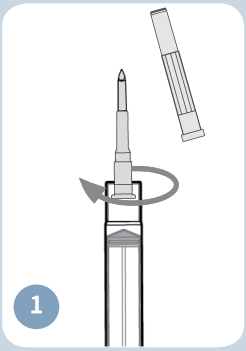


Connect the aspirating needle to the syringe and remove the protective cap

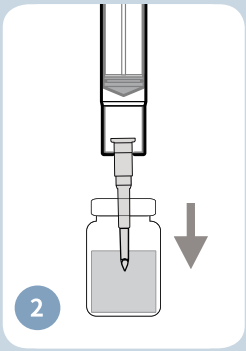


Insert the aspirating needle into the vial

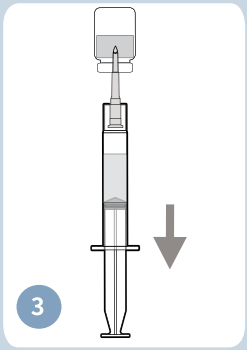


Invert the vial and aspirate the necessary dose of drug liquid for treating patients (including the volume of dead space)

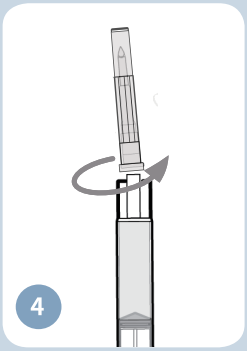


Put on the protective cap and screw off the aspirating needle

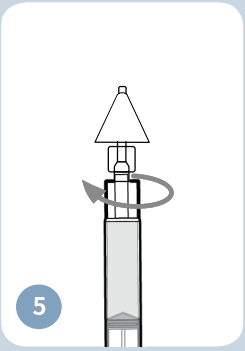


Connect the nasal nozzle to the syringe

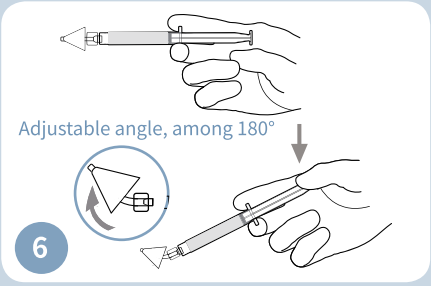

Adjust the angle of nasal nozzle according to the body position of patient

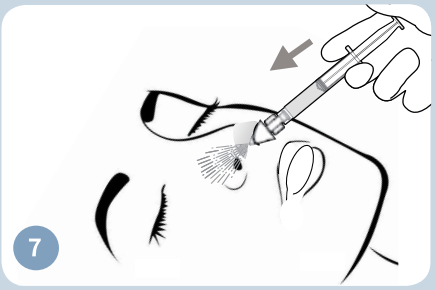

Quickly push the rod, spray half of the drug into the nasal cavity, and complete the nasal administration on one side

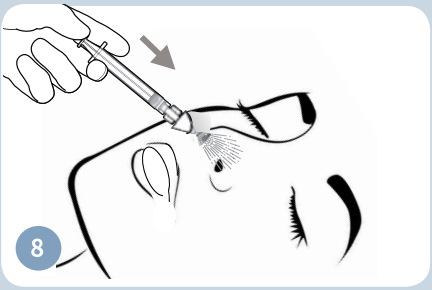

Quickly push the rod, spray the remaining drug into the nasal cavity, and complete the nasal administration on the other side

Specifications

Product	Spray Nozzle Specifications	Injector Syringe Specifications	Puncture Needle type
Single Use Dispensing Nebulizer	(Optional) 5cm/12cm/ 23cm/25cm/30cm	(Optional) 1mL 3mL	(Optional) Side hole needle Bevel needle



Nextech-X

Hangzhou Qiantang Longyue Biotechnology Co., Ltd.
Rooms 302 and 104, Building 12, Block 1, No.619 Wangmei Road, Linping Street, Linping District, Hangzhou City
Website: www.nextech-x.com/en
Email: Sales1@Nextech-x.com
Tel: 0571-89150121, 18017582758